

# Invacare® Bariatric Bed

## Specifications

(all specifications are approximate)

BAR600 Bariatric bed

### Dimensions

Bed height*	18" min. - 26" max.
Sleep surface length	80"
Sleep surface width	42"
Overall length	88"
Overall width	42"
Head-end height	31.5" min. - 39.5" max.
Foot-end height	29" min. - 37" max.
Split spring to ends	43" head - 44.5" foot
	Carton size
	Shipping weight

BAR5000	Head spring	47.25" l x 10" w x 43.75" d	100 lb.
BAR5490	Foot spring	43.5" l x 10" w x 43.75" d	131 lb.
BAR5401	Bed ends	43" l x 9.5" w x 25.5" d	65 lb.

### Limited warranty

All welds	lifetime
Frame/structure	5 years
Mechanical/electrical	1 year

Client weight capacity 600 lb.

\*Measurement from floor to bed deck.




## Model no. BAR600

The Invacare® BAR600 is a heavy-duty bed frame designed for bariatric individuals and capable of supporting up to 600 pounds. Extra bracing has been integrated into the home care bed frame, along with a wider surface and truss assembly, to provide maximum support. The client or caregiver is able to easily operate the bed using the hand pendant, which changes the positioning of the head and foot sections as well as the bed height. Even though this is a substantial bed, setup within the home is easy.

### FEATURES

- 600 lb. maximum client weight capacity
- Full electric functionality
- Heavy-duty bed sleep surface is 17% larger to ensure comfortable positioning
- Split-spring design with removable bed ends for ease in setup
- No special tools necessary for in-home setup
- Easy-release motor mounts allow quick motor removal
- Robotically welded, steel frame construction
- Pendant allows for dynamic client positioning
- Quiet, smooth operation
- Emergency crank manually functions to lower any raised section of bed
- UL73 approved
- Use with Invacare bariatric mattress (BARMATT) and bed rails (6640)
- Designed to meet KO549

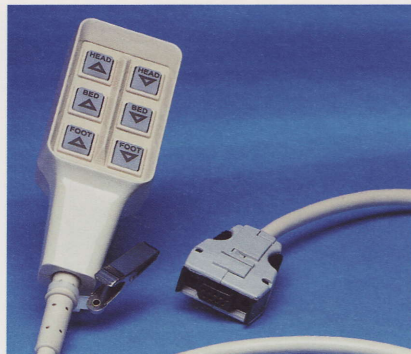


Hand pendant.



BAR600 with head raised.

## Additional Product Suggestions

Additional product suggestions are for recommendation only. Please refer to the individual product specification to verify weight capacity.

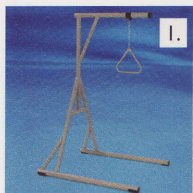

1.



2.



3.

1. Bariatric Trapeze #BARTRAP, page 26
2. Reliant 600 Heavy-Duty Lift #RPA600-I, page 27
3. Bedpad Reusable Waterproof 34 x 52 in. #ISG13452, call Invacare Supply Group at 800-225-4792

Invacare Customer Service (800) 333-6900 [www.invacare.com](http://www.invacare.com) Invacare Bariatric Price List Form no. 99-277